


Technical Data

ZetaView® Evolution mono

PARTICLE
MATRIX

analytik

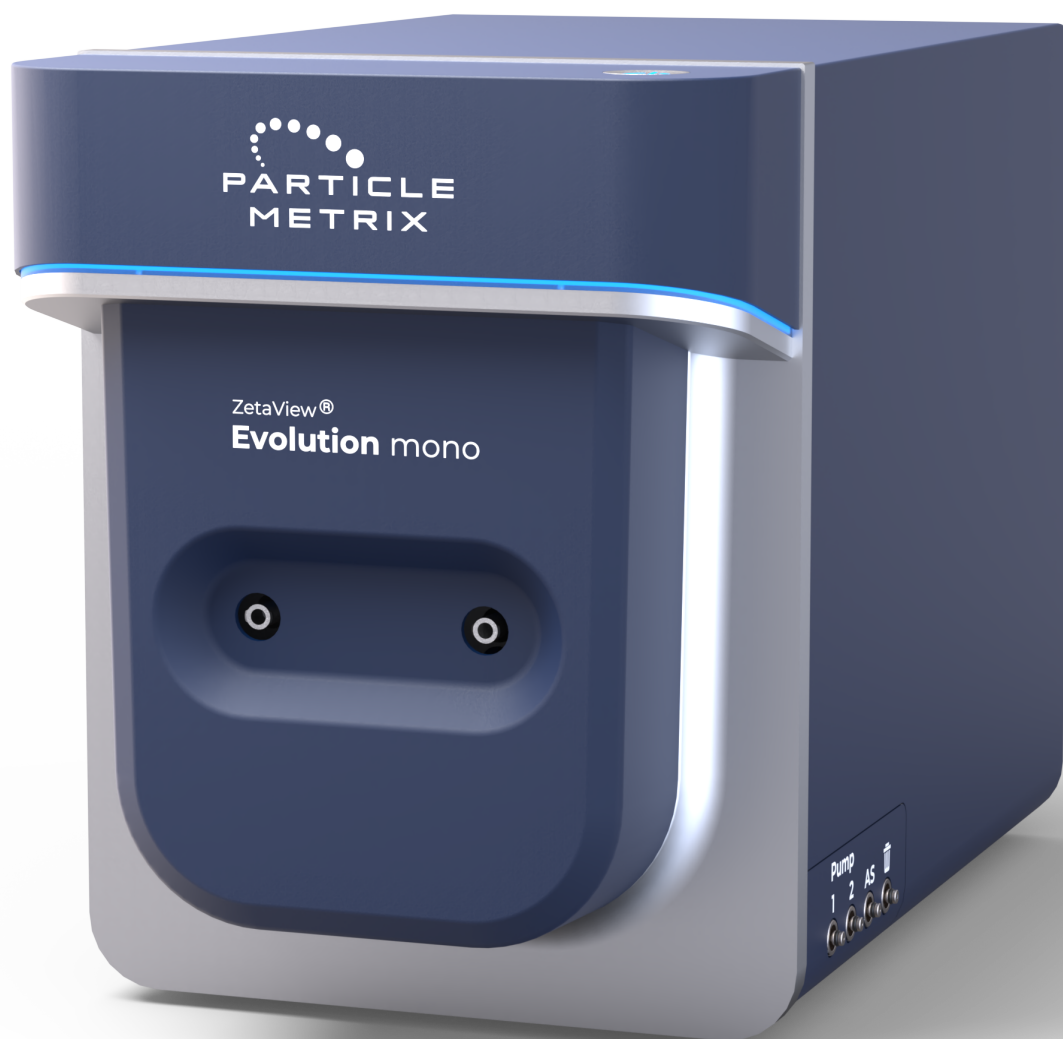

Size



Fluorescence



Colocalization



Zeta
Potential



Sub-
populations



Concentration



Technical Data

ZetaView® Evolution mono



Measurement principle in scatter and fluorescence

Size	<ul style="list-style-type: none">Motorized scanning Nanoparticle Scanning Analysis (NTA) for measurement of up to 39 subvolumes per sample (129nl measurement volume)
Concentration	<ul style="list-style-type: none">Concentration Scanning Technology or NTA (220nl or 129nl measurement volume)
Zeta potential*	<ul style="list-style-type: none">Video based electrophoretic mobility tracking

Data management

Software	<ul style="list-style-type: none">ZetaSphere control software featuring measurement of size, concentration and zeta potential* in scatter and fluorescence
Quality control	<ul style="list-style-type: none">Integrated instrument performance checkOutlayer control by automatic Grubbs statistical analysis of measurement dataDatabase event logging for data integrityLive monitoring of particle size and concentration, temperature, scattering intensity, conductivity, particle drift and signal to noise ratioPredefined measurement settings for several applications which are fully customizable
Data output	<ul style="list-style-type: none">PDF multiparamter sample reportingCSVPNGFCS

Physical characteristics

Dimensions	<ul style="list-style-type: none">W x D x H: 20cm x 30cm x 25cm (excluding computer)
Weight	<ul style="list-style-type: none">13.5kg (excluding computer)
Shipment conditions	<ul style="list-style-type: none">on pallet
Electrical supply	<ul style="list-style-type: none">90V - 240V; 47 - 63Hz; 50VA
Power consumption	<ul style="list-style-type: none">max 30W

*When ordered with zeta potential option

Technical Data

ZetaView® Evolution mono



Hardware

Instrument

- 90° laser scattering video microscope with x10 magnification for maximized sample volume and highest statistics
- One software controlled laser for use in scatter and fluorescence
- Software controlled 12 position fluorescence filter wheel for scatter and 11 fluorescence channels
- Two software controlled pumps for liquid transport and sample dosing
- Designed for automated sample loading
- Automated alignment and focussing of laser and microscope
- External temperature range: 5°C to 45°C
- Sample temperature control via peltier element from RTP -5°C to 55°C with automated due point sensing

Camera

- High sensitive CMOS camera with 1280x960 pixels
- Variable frame rate from 2Hz to 60Hz for optimal resolution and fast acquisition

Lasers

- Available laser options: 405nm (130mW); 488nm (40mW); 520nm (80mW) or 640nm (135mW)
- Pulse duration each laser: 0.1ms up to continuous

Filters

- Software controlled automated 12 position filter wheel equipped with one fluorescence emission long pass filter at 430nm (in case of 405nm laser), 500nm (in case of 488nm laser), 550nm (in case of 520nm laser), and 670nm (in case of 640nm), cut-off
- customized emission filters available on request

Measurement cell

- Quartz class cuvette for low protein binding
- Tool free access for quick and simple cleaning process

Computer

- i5 Asus® NUC Mini PC (i7 optional)
- 1TB SSD hard drive
- Windows 11 Professional
- Keyboard and mouse

Monitor

- LED screen

*When ordered with zeta potential option

Technical Data

ZetaView® Evolution mono



Measurement specifications

Size	<ul style="list-style-type: none">• 10 – 1000nm (dependent on sample and laser)• Accuracy: $\pm 5\%$ (for 100nm polystyrene latex)
Concentration	<ul style="list-style-type: none">• Concentration range: 10^5 – 10^9 particles/ml• Accuracy: $\pm 5\%$ (for 100nm polystyrene latex)
Fluorescence	<ul style="list-style-type: none">• Concentration range: 10^5 – 10^9 particles/ml• Sensitivity level: < 20AF488 molecules
Zeta potential*	<ul style="list-style-type: none">• Working range: -500mV - +500mV• Concentration range: 10^6 – 10^{10} particles/ml• Conductivity range: $3\mu\text{S/cm}$ – 15mS/cm• Reproducibility: $\pm 2\text{mV}$ (zeta potential standard)
General	<ul style="list-style-type: none">• Minimum sample quantity: 500μl of sample at 10^5 particles/ml• pH range: 1 - 13
Reference material	<ul style="list-style-type: none">• Nominal 100nm size and concentration reference suspension• Four nominal 100nm and 200nm reference suspensions for fluorescence• Nominal -50mV reference suspension for zeta potential*

*When ordered with zeta potential option

Manufactured by Particle Metrix, distributed in the UK and Ireland by **analytik**.